

Designed and manufactured in Australia by Thermoline

**ABN:** 80 000 859 129

**Head Office:** 10-12 Ross Place Wetherill Park NSW 2164 Australia **Phone:** +61 2 9604 3911

Email: hello@thermoline.com.au

Web: www.thermoline.com.au

# Proudly Family-Owned and Operated Today, we remain a proud second-generation family business





#### Designed and Manufactured by Thermoline

The Thermoline Laboratory Drying Oven models are designed and manufactured to provide safe and precise temperature conditions for the items stored inside.

Thermoline's Australian Made drying ovens are designed to remove moisture from products and samples. Adjustable air inlet and exhaust vents supply fresh air into the chamber while removing moist air from the chamber.

Typical applications include moisture removal from soil, paper, timber, plaster, fabric, coal and other minerals.

#### Models available:

Thermoline have six drying ovens available, all designed and manufactured in Australia. These models are economical and very easy to use making them an ideal choice for all processes where products and samples require drying.






#### • Digital Controller:

The Omron E5CC-T uses an automated programmable microprocessor to control the temperature within +/-0.1°C. The bright LED display has a range of functions including high alarm, temperature calibration and programmable multi-setpoint functions.

#### Fan Forced Air Circulation:

Fan-forced air circulation is a process of circulating air using a fan to distribute it evenly throughout the space. Thermoline's horizontal airflow blows dry air across each shelf for effective drying anywhere inside the oven. Heated air is blown down the rear of the oven where it is then forced over the shelf spaces via a perforated rear wall.

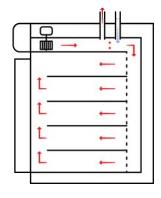



#### Stainless Steel Interior:

Stainless steel is durable, easy to clean and corrosionresistant, making it the perfect choice for high heat tasks. 316 marine-grade stainless steel is used here to ensure the most corrosive-resistant liner possible.


#### Australian Made:

This range is designed and manufactured by Thermoline. We proudly promote and support the Australian Made logo.





#### Designed and Manufactured by **Thermoline**

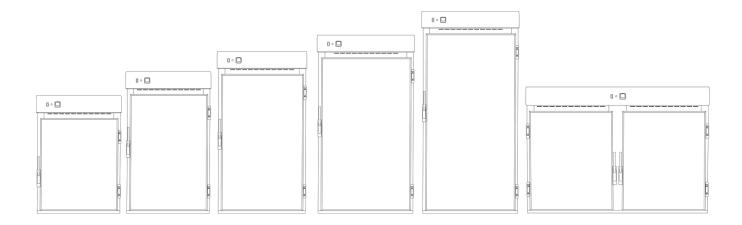

#### **Dimensions**


External	TD-80F	TD-150F	TD-250F	TD-500F	TD-700F	TD-630F
WxDxH (mm)	520x560x800	630x660x1010	630x660x1360	780x810x1680	780x810x1980	1450x860x1100

Note: Inlet and Exhaust ports extend the overall height of the ovens by approximately 40mm.

#### Internal

WxDxH (mm)	390x400x400	500x510x600	500x510x950	650x650x1200	650x650x1500	1300x690x650
------------	-------------	-------------	-------------	--------------	--------------	--------------

Clearance	TD-80F	TD-150F	TD-250F	TD-500F	TD-700F	TD-630F
Front (mm)	520	630		780		725
Back (mm)	100					
Sides (mm)	100					

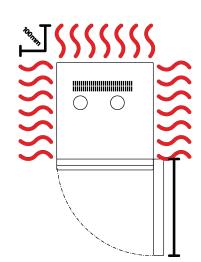



See our Volatile Materials Oven Range. The TD-AS1681.

#### Designed and Manufactured by **Thermoline**

Technical Specification	TD-80F	TD-150F	TD-250F	TD-500F	TD-700F	TD-630F	
Temperature Range	Ambient +10°C to 200°C						
Temperature Control Stability	+/- 0.1°C						
Temperature Uniformity		+/-2°C @ 100°C		+/-4°C @ 100°C			
Heater Power		1500 watts		3000 watts			
Electrical		7A/230V		13A/230V (requires 15A supply)			
Nominal Capacity	80L	150L	250L	500L	700L	630L	
Porthole Diameter	13mm						
Weight	60kg	100kg	110kg	170kg	190kg	190kg	



## Expert Support Team

Our R&D, manufacturing, support, and sales teams are all in-house allowing for efficient collaboration


Designed and Manufactured by **Thermoline** 

Features	TD-80F	TD-150F	TD-250F	TD-500F	TD-700F	TD-630F
Shelves (max @100mm spacing)	2 (max 3)	3 (max 4)	4 (max 6)	5 (max 9)	6 (max 11)	3 levels (max 4)
Lockable Castors	Optional	/	/	/	/	/
Internal Fan	/	/	/	1	/	✓
Omron E5CC-T	/	/	/	/	/	<b>✓</b>
Solid Door	/	/	/	/	/	<b>✓</b>
Fibreglass Insulation	/	/	/	1	/	✓
Safety						
Over Temperature Safety	/	/	/	•	/	/
<b>Over Current Protection</b>	/	/	/	/	/	~

#### **Options**

BMS Plug	No volt contact closure plug and socket connection to a Building Management System
Additional Shelves	Additional Stainless Steel shelves to suit
Heavy Duty Shelves	Update to allow for higher loads on shelves (limits total number of shelves possible for use)
Door Locks	Key lockable door locks
Customisable Port Hole	Add additional stainless steel 13mm port holes or choose 50mm port holes
Programmable Controller (MSP)	(E5CC-MSP) Upgraded Omron Controller features 256 ramp and/or dwell segments
Glassware Drying	See TGD-range for oven models specific to drying of laboratory glassware
Natural Convection	See TO-range for oven models without fan forced circulation



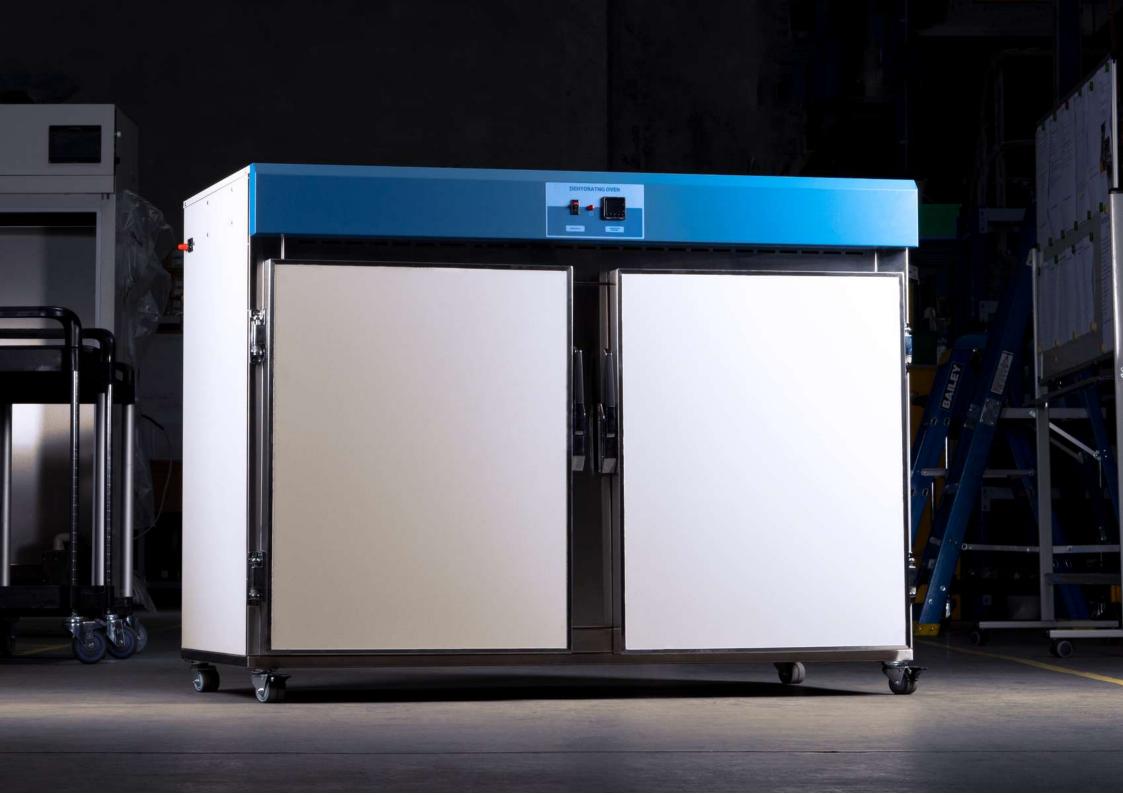

## We are proudly Australian owned

We will continue to invest in Australian manufacturing.

